

# Preparations for Use




Dry Volume	Recommended BMA Volume
5cc	2.5cc
10cc	5cc
15cc	7.5cc

1. Remove the putty from the sterile packaging.

**Note:** Before hydration, putty is formed in a flat dry block.

2. The dry putty should be hydrated with Bone Marrow Aspirate (BMA) using approximately  $\frac{1}{2}$  of the dry volume to hydrate the product. (See Chart).

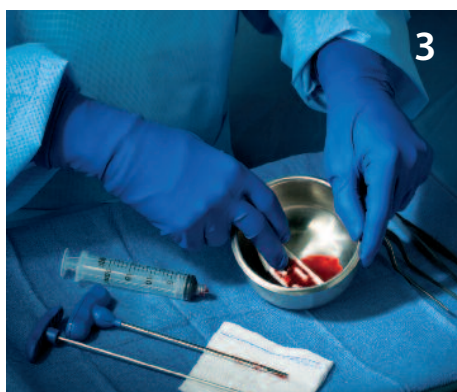

3. Dry putty absorbing BMA.

4. Fully Hydrated putty.



5. Once fully hydrated and mixed, the product will soften and become moldable into the desired shape for subsequent placement into the surgical site or operative defect. Additional fluid can be added to achieve desired consistency.

\* Refer to package insert