




RIVAL™

○ REDUCE

□ VIEW

▶ BITE

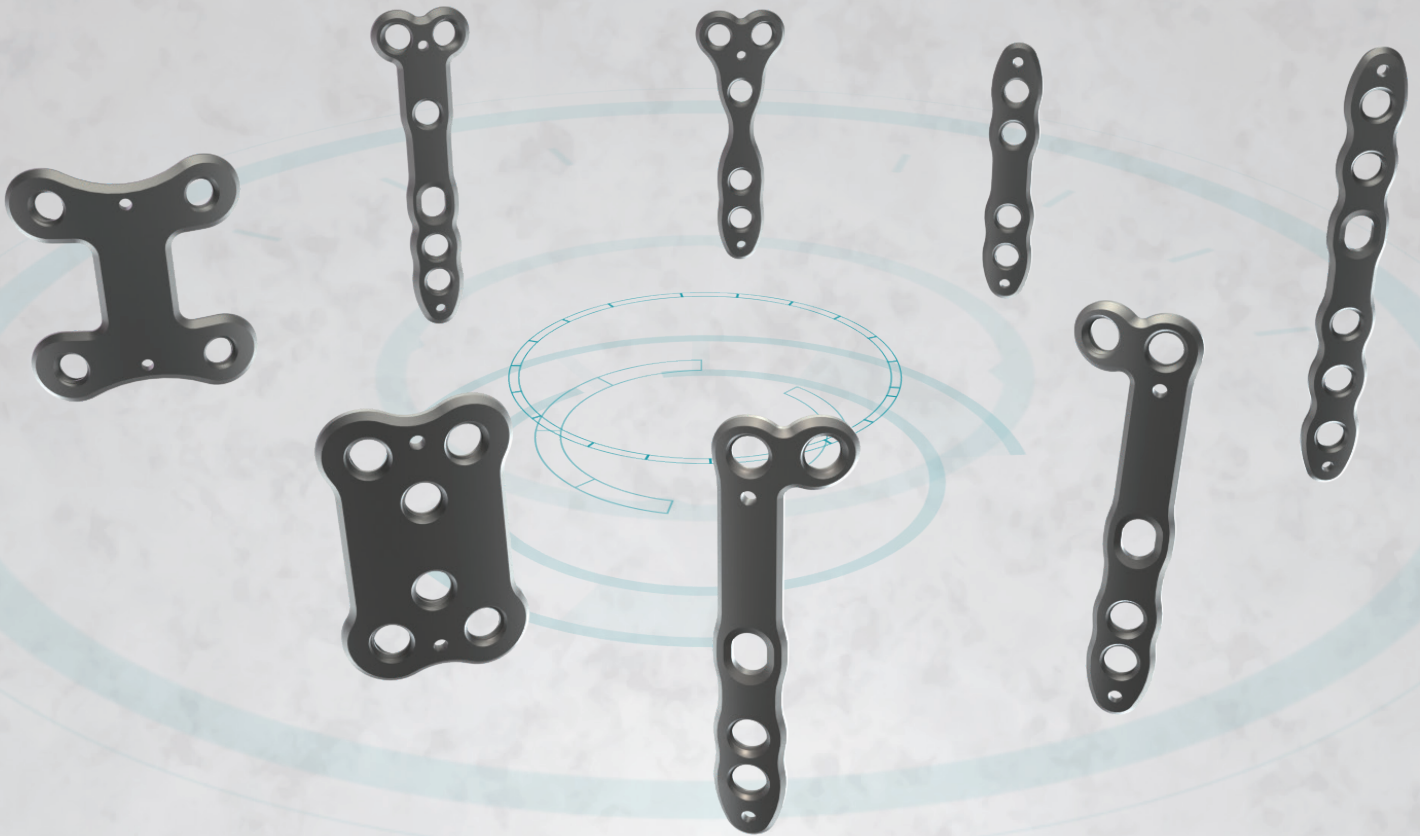

○ REDUCE


Fracture Plating System

Wide range of
implants and sizes

Optimal implant selection


Designed to
reduce soft tissue irritation

Designed to
reduce screw prominence

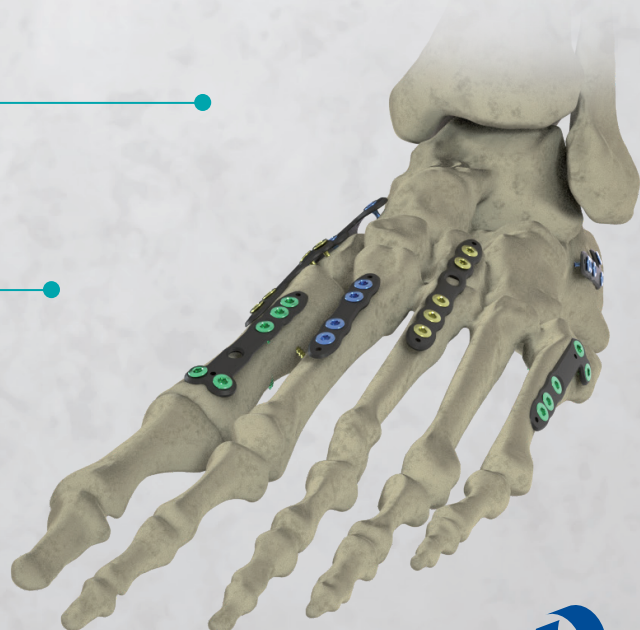

Designed to be used with 2.4mm, 2.7mm and 3.2mm locking and non-locking screws

Locking screws

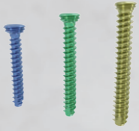

Non-locking screws




Guide wire holes

For easy implant placement during surgery



Double lead thread screw design

For fast and controlled insertion

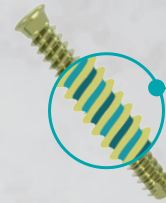

Oblong compression hole

Intraoperative compression up to 1.1mm with the corresponding non-locking screws

Type II anodised titanium for **Increased fatigue resistance**



Low profile

1.3mm plate thickness





Screws lock flush into plates reducing screw prominence


Universal instrumentation

Designed


- to be used across all RIVAL™ product lines
- to help with inventory management
- to provide a comprehensive range of instruments

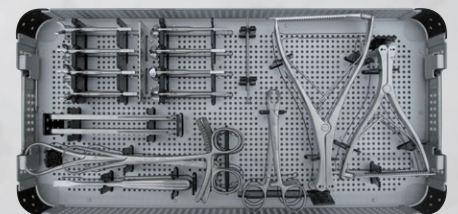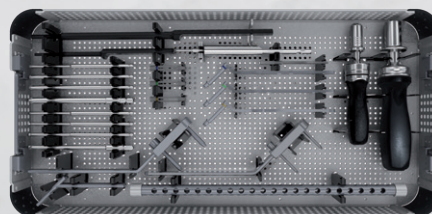

Distributed by:

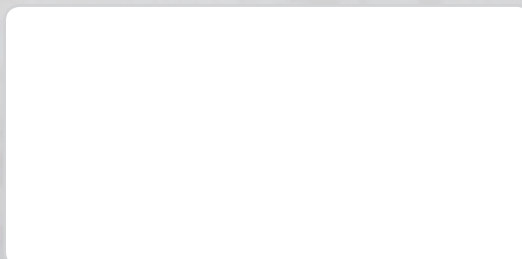

Manufactured by:
 ORTHOFIX Srl
 Via Delle Nazioni 9, 37012 Bussolengo (Verona), Italy
 Telephone +39 045 6719000, Fax +39 045 6719380

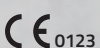

www.orthofix.com

RIRE-1904-PL-E0 A 10/19

